
MV11 Portable Q-Switched Nd:YAG Laser Pigment Lesions & Tattoo Removal Beauty Salon Machine

HONKON
Science Creates Beauty



HONKON
Science Creates Beauty



HONKON
Science Creates Beauty



HONKON
Science Creates Beauty



HONKON
Science Creates Beauty



Powerful Portable Q-Switched Nd:YAG Laser Pigment Lesions & Tattoo Removal Beauty Salon Machine

Powerful portable 1064nm & 532nm Q-switched Nd:YAG laser

Pigment lesions & tattoo removal

MV11



◆ Safe and efficient

◆ High Cost-effective

◆ Easy to operate

◆ Significant effect

Treatment Therapy


The Q-switched Nd:YAG laser radiates the light with a specific wavelength in high peak energy pulses, hence the light penetrates into the tissue only for a nanosecond. The light is absorbed by the pigmentation and results in an instantaneous blast, that is Light blasting principle. The pigmentation particles are shattered into fragments, some can be bounced out of the skin and others can be split into tiny particles that can be engulfed by phagocytes and then eliminated by the lymphatic system.

MV11 TREATMENT PRINCIPLE

Treatment principles

The Q-switched Nd:YAG laser delivers light of specific wavelengths in ultra-short (nanosecond). The high peak energy pulses are absorbed by the pigments and results in an instantaneous blast, that is light blasting principle. The pigments and results are shattered into fragments, some can be bounced out of the skin and others can be split into tiny particles that can be engulfed by phagocytes and then eliminated by the lymphatic system.

Process of laser removing pigment

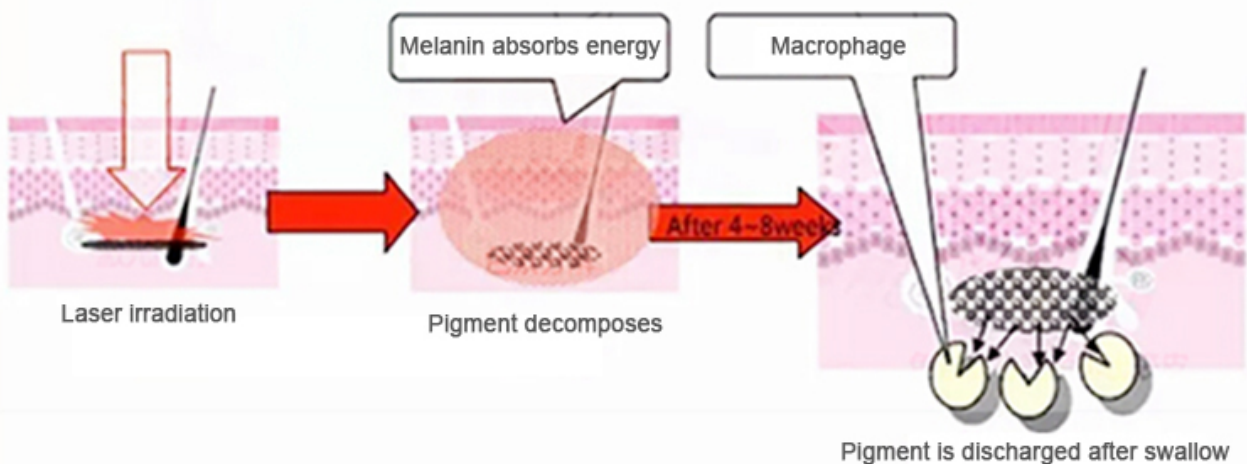

Application

1. Tattoo (black, brown, red, blue ,etc.) Tattoo eye line,tattoo eyebrow
2. Pigmented lesions : Sunspot ,Age spots, Nevus,Freckle.
3. Laser facial : Skin rejuvenation , fine wrinkles and expression lines removal, pore firming,lighten spots,oil control

Tattoo Removal ND: YAG Q-Switched Laser

Application

OWNER DEDICATED AIR PURIFICATION
www.foxconn.com



Exogenous pigmented skin problems

Tattoo, Eyebrow tattoo, Eyeliner tattoo, Lip line and traumatic pigmentation

Endogenous pigmented skin problems

Freckle, Age spot, Black mole etc.



With the toner treatment head MV11 can easily carry out the skin rejuvenation, to achieve rejuvenation, shrinking pores, fading the stain, skin whitening, blackhead removal, facial oil control.

Tattoo Removal ND: YAG Q-Switched Laser

Structure display



Intelligent laser handles

1064nm and 532nm treatment head easily replaced

此设备为医疗美容设备，使用时请严格按照说明书操作。本产品为二类医疗器械，使用时请严格按照说明书操作。本产品为二类医疗器械，使用时请严格按照说明书操作。本产品为二类医疗器械，使用时请严格按照说明书操作。



1064nm Treatment head



532nm Treatment head



Carbon powder(Skin rejuvenation) treatment head



Features & Advantages

1. MV11 is with keyboard screen
2. Patent technology of single lamp and dual rods doubles energy output(400mj).
3. Shorter pulse width(6-8ns)than ordinary with carbon skin rejuvenation function.
4. Potable design and easy to remove,the cost of use and maintenance is lower.
5. Spot adjustor is optional,high cost-performance.
6. Better results for dermal pigmentation

Treatment Advantage



More comfortable

The laser whose width less than 8ns shortens the working time on tissue and reduces the pain.

Easy to carry

Tiny and portable, easy for transport and recycle therapy

Faster

Unique single lamp double tubes patented technology, single pulse energy with over 400mj is twice than normal equipment, treatment time is reduced by twice.

Specification

Laser-Type	MV11	Working frequency	1-5Hz
Wavelength	1064nm/532nm	Cooling mode	Water-cooling + air-cooling
Single pulse energy	≤400mj@1064nm ≤200mj@532nm	Spot diameter	1~6mm Continuously adjustable
Double pulse energy	≤800mj@1064nm ≤400mj@532nm	Package dimension	535*532*745mm ³
Triple pulse energy	≤1200mj@1064nm ≤600mj@532nm	Gross Weight	20kg
Pulse width	6-8ns	Insurance specifications	Max5A
Power supply	AC110 - 220V ±10%, 50Hz - 60Hz	Input power	780W

10:11

4G



She want to buy

Yes this is ours best selling I will send you some documents



So please quote me your best price first

Rmb?



Or usd



Usd please

1064QCCL.pdf
5M



Instruction Maual of
1064QCCL.pdf
631K



I have 1 Spa ,too. I used your machine a longtime



Very good quality



Product link : <https://www.hongkongglobal.com/mv11.html>